Reed Biotech

VeriQuant™

All Components Included
Get Started Immediately.

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
FULL-STACK SELF-DEVELOPED SYSTEM EMPOWERING COMPREHENSIVE PERFORMANCE ENHANCEMENT

(a) Two monoclonal antibodies against distinct epitopes are conjugated with oligo probes, then form a sandwich complex with the antigen.

(b) Following antigen recognition by the paired antibodies, the conjugated oligonucleotides are spatially proximate, and generate an amplifiable dsDNA template under the enzymatic action.

(c) qPCR amplification is then performed. CT values correlate inversely with protein concentration, and sample target concentrations are calculated via standard curve.

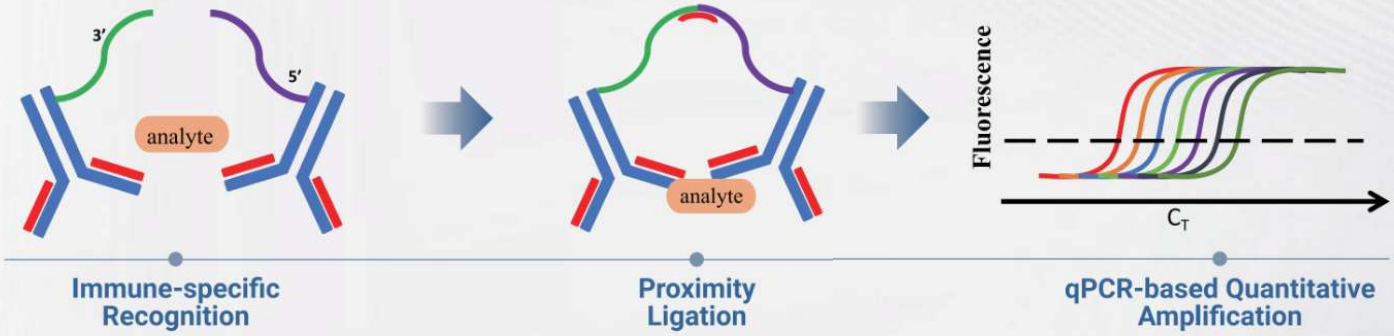Reed's full-stack self-developed polymerase Quanti Mixtureand upgraded Synthesis Enzymesbuffer system can easily achieve fg/mL-level protein quantification; Greatly simplify the process, enhance efficiency and consistency. Integrate the "connection" and "extension" reactions in the same reaction tube, and complete all chemical reactions with just one incubation without the additional pipetting.